

TheraFit Premium Soft Goods

Shoulder Harness - Vest - Chest Strap - **4 Point Hip Belt** - Hip Belt - Rear Pull Hip Belt - Center Pull Hip Belt



Custom Seating 2-25

TheraFit Soft Goods

Trays 77-98

Tray Attaching 99-103

Sip-N-Puff / Joystick 104-107

Patient Transfer 108-111

Therapy 112-115

Independent Living 116-121

TheraFit 4 Point Hip Belts

- Choose 1" or 1.5" straps, adjustable or fixed
- Various pad sizes available
- Metal push button buckle or plastic side release buckle
- Hardware to attach to chair included

TheraFit 4 Point Hip Belt lengths shown here are taken from end to end, fully extended, without attaching hardware.

Belt w/1" Straps, 4.25"L x 1.75"H Pads

w/Adjustable Straps & Metal Buckle, 44"	.60125
w/Adjustable Straps & Plastic Buckle, 44"	.60135
w/Fixed Straps & Metal Buckle, 48"	.60130
w/Fixed Straps & Plastic Buckle, 48"	.60140

1" STRAPS

Belt w/1" Straps, 5.25"L x 1.75"H Pads

w/Adjustable Straps & Metal Buckle, 46"	.60126
w/Adjustable Straps & Plastic Buckle, 46"	.60136
w/Fixed Straps & Metal Buckle, 50"	.60131
w/Fixed Straps & Plastic Buckle, 50"	.60141

Belt w/1.5" Straps, 5.25"L x 2.25"H Pads

w/Adjustable Straps & Metal Buckle, 51"	.60127
w/Adjustable Straps & Plastic Buckle, 51"	.60137
w/Fixed Straps & Metal Buckle, 50"	.60132
w/Fixed Straps & Plastic Buckle, 50"	.60142

1.5" STRAPS

Belt w/1.5" Straps, 7.25"L x 2.25"H Pads

w/Adjustable Straps & Metal Buckle, 55"	.60128
w/Adjustable Straps & Plastic Buckle, 55"	.60138
w/Fixed Straps & Metal Buckle, 54"	.60133
w/Fixed Straps & Plastic Buckle, 54"	.60143



This example has 1" Straps that are adjustable

Metal Push Button Buckle



This example has 1.5" Straps that are fixed

Plastic Side Release Buckle


HIP BELTS & CHEST STRAPS

*Parts included vary depending on which belt you have.
You may not have all parts shown here.*


METAL ATTACHING CLIPS

Depending on width of strap, thread strap as shown.

Metal Attaching Clip & 1" Or 1.5" Strap

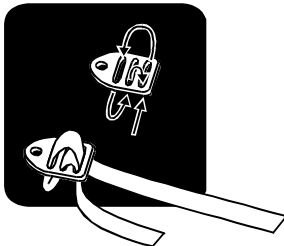


Metal Attaching Clip & 2" Strap

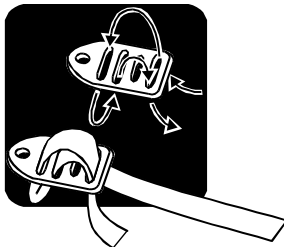


CAM LOCK BUCKLES

Thread strap through buckle and push top down to hold strap in place.

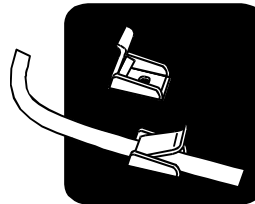


SLIDING DOUBLE-D CLIPS

Use Sliding Double-D Clips to loop straps around chair frame.

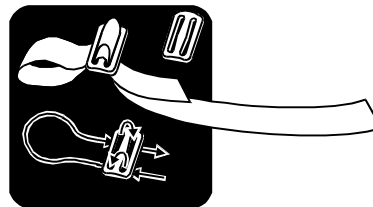


Use Sliding Double-D Clips with Metal Attaching Clips.

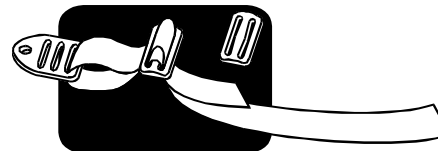

HINGED COLLARS

Attach Hinged Collars to chair frame. Attach Metal Attaching Clips or Cam Lock Buckles to Hinged Collars.

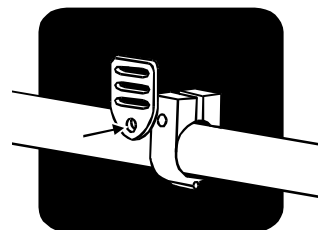

08756