



HydraTherm & HydraHeat



The next generation of moist heat therapy....

A revolutionary composite heating unit that breaks the boundaries of traditional moist heat therapy. Introducing the Performa HydraTherm, the first ever rust, mold, bacteria and mildew-proof unit, that allows for a clinical environment that is:

✓ CLEANER

✓ SAFER

✓ MODERN

...addresses your largest frustrations.

Infection Control Issue

Industry mandated cleaning and temperature logging

- ✓ High grade polymer prevents cross contamination by eliminating bacterial growth
- ✓ Digital read out for precise temperature monitoring
- ✓ New design significantly reduces the number of annual cleanings (2 versus 24)

Safety Issue

Patient burns and facility fire hazards

- ✓ Digital thermostat enables clinician to control treatment temperature
- ✓ Low water sensors activate automatic shut off to reduce fire hazards
- ✓ 24V versus 110V decreases electrical safety threats

Dated Technology

Lack of innovation

- ✓ Automatic pump drains water at the flip of a switch
- ✓ Divider racking system provides easy access to hot packs
- ✓ Lower voltage increases energy efficiency





Divider Racking System

Easier access to hot packs eliminates the need for tongs

Fold-Away Shelf

(Deluxe Model)

Fold-Away Shelf is easily positioned to create additional work area or lowered for space utilization



Automatic drain pump

(Deluxe Model)

Easily remove water with the flip of a switch

Ordering Details

REF	Description
081686088	Hydratherm Standard - Divider Rack System Includes 3 Standard, 3 Universal, 3 Cervical, 3 Oversize Performa HydraHeat Packs
081686096	Hydratherm Standard - Traditional Rack System Includes 3 Standard, 3 Universal, 3 Cervical, 3 Oversize Performa HydraHeat Packs
081686120	Hydratherm Deluxe - Divider Rack System Includes 3 standard, 3 Universal, 3 cervical, 3 oversize Performa HydraHeat Packs, Fold-Away Shelf, Automatic Drain Pump
081686138	Hydratherm Deluxe - Traditional Rack System Includes 3 Standard, 3 Universal, 3 Cervical, 3 Oversize Performa HydraHeat Packs, Fold-Away Shelf, Automatic Drain Pump

Technical Specifications

Power	110/120 VAC - 220/240 VAC 50/60 Hz
Power consumption	HT-R12-S/HT-R12-D 1000 Watts
Adjustable temperature range	120°F (49°C) 160°F (71°C)
Safety thermal cut-out temperature	175°F +/- 5°F (79°C +/- 3°C)
Digital temperature control accuracy	+/- 5%
Weight (dry)	
Standard model	70 lbs (31kg)
Deluxe model	75 lbs (34kg)
Water capacity (without heat packs)	18 Gal (68L)
Dimensions	30"W x 20"D x 33" H (76cm x 50cm x 84cm)

Performa HydraHeat Packs



Innovative technology that delivers therapeutic moist heat at safer temperatures without the worry of bursting. This latex free formulation is compatible with all types of heating units.

Safety Solution

Lower required temperature

- ✓ Reduces the risk of patient burns
- ✓ Eliminates the need for terry covers or towels

Infection Control Solution

Non-porous, bacteria resistant, polymer surface

- ✓ Allows for easy sanitation
- ✓ Prevents unit and patient cross contamination




Quality Solution

Patented heat absorbing gel formula with SureSeam technology


- ✓ Protects against bursting and leaking
- ✓ Provides even heat distribution with greater surface contact

Ordering Details


REF	Description
081686179	Standard 17"x 11" (43 x 28cm)

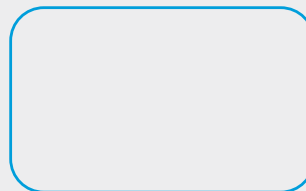


081686161	Neck Contour/Cervical 17"x 11" (43 x 28cm)
-----------	--

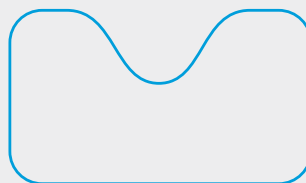

081686153	Universal 19"x 5" (48 x 13cm)
-----------	---

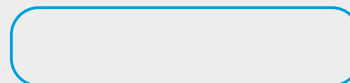

081686146	Oversize 17"x 13" (43 x 33cm)
-----------	---

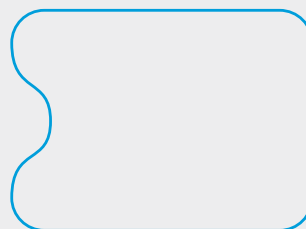

Call Customer Service: 1-800-323-5547
24-Hour Fax: 1-800-547-4333
Customer Service Hours: 7:00am - 7:00pm CST

081687789